

Response™

from  Comfortex

Stretcher Cushions

Comfortex Response™ Stretcher Cushions

A comfortable, well engineered, stretcher cushion goes a long way towards providing your patients with a better healthcare experience.

- * Pressure Reducing Protection
- * Comfortable Dual Layer, High Resiliency Core
- * Replaceable Top Covers (shear reducing Stratos Nylon or Lectrolite Vinyl)
- * Embedded CPR boards



 Comfortex

800-445-4007

*Being a Patient can be stressful enough.
Reduce Pressure ... Call Comfortex*

Response™

from  Comfortex.

Stretcher Cushions

Design Features

Mattress Construction & Comfort by Comfortex, one of the nations leading manufacturers of pressure reducing surfaces.



Replaceable Components allow cushions to return to service in quick order should damage occur. Whether to replace a cut cover or simply to restore appearance, zip off top covers can be replaced for just \$48.

Loose Fitted Covering allows the cushion to conform to the patient's unique shape & weight. It looks a little different, but feels great.

Shear Reducing Stratos™ nylon top covers are healthcare's most durable barrier against fluid penetration. A special Lectrolite anti-static vinyl is also available. (FYI : Stretch & breathable fabric are not used, as they can allow moisture pass thru and cause unseen core contamination)

Embedded CPR Boards are undetected by the patient or by x-rays. A cardiac arrest can occur anywhere, anytime...this cushion Responds.




Design Features, cont.

Never a pillow when you need one?

Fiber filled pillows can be ordered and attached to your stretchers. Flip in place when desired or flip out of the way if the patient needs to lay flat. Webbing straps and snap-clips secure the pillow to the stretchers frame.



Need Affordable Excellence...Comfortex responds.

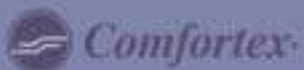

800-445-4007

***Being a Patient can be stressful enough.
Reduce Pressure ... Call Comfortex***


Response™

from  Comfortex.

Stretcher Cushions

Clinical Performance

Whether a short trip from room to radiology, or a pre-op, post-op marathon, the stretcher cushion can have a huge impact on your patient's healthcare experience.

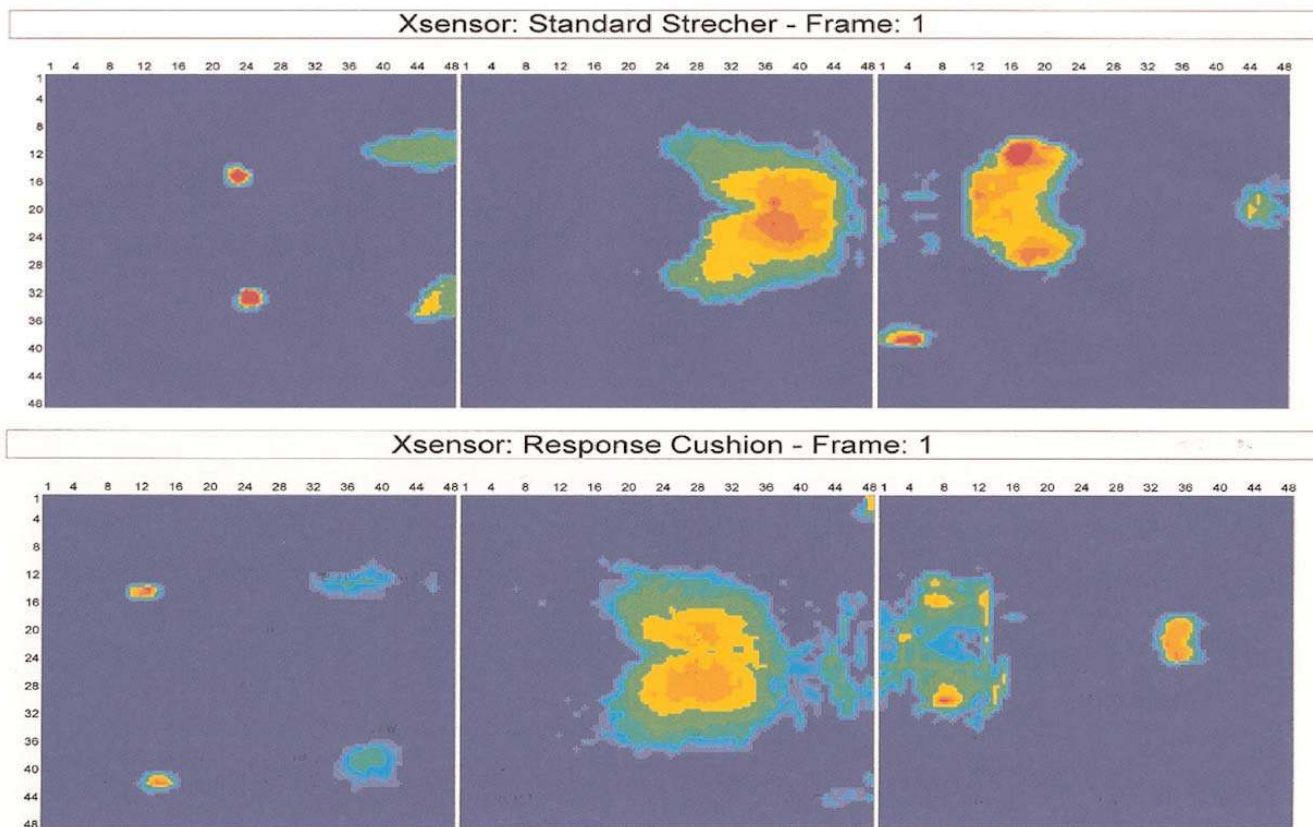

When compared to conventional pads, the Response Cushion better conforms to the patients shape, weight & movements for better pressure redistribution & greater comfort.

We've all sat on a hard, poorly cushioned surface and after just a short time, began to shift & reposition in an effort to find comfort. However, many times patients aren't able to shift themselves.

With our Response Stretcher Cushion we've blended our most durable shear reducing fabric, sewn in a loose, billowy fashion and HR "High Resiliency" core materials to produce a cushion that supports the body and comfortably cushions sensitive areas.

Computerized Pressure Imaging illustrates how the Response Cushion reduces interface pressure at the areas most problematic. **Even when on the cushion short-term, these reductions in localized pressure will help to prevent the discomfort and anxiety patients can experience.**

Clinical Performance, cont.

Memory Foam...Don't get the wrong Impression

Comfortex vs. Visco - Memory Foam

The German Institute of Healthcare Research studied the effect two of Europe's leading pressure reducing mattresses would have on patients who were suffering from chronic pain.

The results were as follows:

All but 2 of the 80 patients reported a change in their level of pain & discomfort when moved from the conventional mattress. The number of patients who awoke as a result of pain was reduced by 21% on the Visco – Memory Foam mattress and reduced by a much greater 52% when on the Comfortex mattress.



The need for pain medications was not reduced for patients on the Memory Foam Mattress. However, pain medication dosages were reduced by 24% during night time hours and by 16.5% during daytime hours for patients on the Comfortex surface.

For better healthcare experiences...Comfortex responds

"Some patients are on these stretcher cushions for many hours and we are able to reduce break down substantially. We won't buy anything else."

Materials Manager - Omaha, NE



800-445-4007

*Being a Patient can be stressful enough.
Reduce Pressure ... Call Comfortex*

Response™

from  Comfortex

Stretcher Cushions

CPR Support

Cardiac Arrest can occur in an elevator as well as in the ER. The Response cushion has x-ray permeable boards embedded beneath the patient's chest area.


Positioned beneath the chest area, between the soft, top layer and supportive base, the CPR boards are virtually imperceptible—unless direct pressure is applied. We hope this feature is rarely utilized, however if need, it's appreciated. (2 boards are in place within non head-end specific cushions)



When seconds count... this cushion Responds.

"We love the fact that the CPR boards are already in the cushion. This saves us a lot of time in an already stressful situation."

ER Supervisor - Springfield, MO

 Comfortex

800-445-4007

*Being a Patient can be stressful enough.
Reduce Pressure ... Call Comfortex*

Response™

from  Comfortex

Stretcher Cushions

Infection Control

Prevent Vapor Flow

Stretch or Breathable fabrics may be waterproof but have High MVT ratings and can permit evaporating fluids to pass thru which can then accumulate in and contaminate the cushion's core.

Comfortex uses only non-stretch fabrics with low MVT ratings to prevent both moisture and moisture vapor pass thru. (see test data : Comtec Laboratories, 2007)



*stretch fabric

Fabric Comparison


Comfortex Stratos Nylon 0.5 to 0.8 Grams/hr/m²

prevents Moisture Vapor Transmission

Waterproof Stretch Fabric 9.4 to 16.0 Grams/hr/m²

permits Moisture Vapor Transmission

For maximum protection...Comfortex responds.

 Comfortex

800-445-4007

*Being a Patient can be stressful enough.
Reduce Pressure ... Call Comfortex*

Response™

from  Comfortex

Stretcher Cushions

Pricing

For a low price and great value...Comfortex responds.


Quantity	1-5	6-11	12+	Contract Pricing
Price	\$258.00	\$250.00	\$240.00	Call Comfortex
Pillows	\$48.00	\$46.00	\$43.00	



Made to Fit

- * For maximum comfort & pressure reduction, all Response cushions are dual layered and 4" thick. (3" available upon request)
- * All Response cushions include embedded X-ray permeable C.P.R. boards and velcro attachment strips.
- * If multiple cushions of various sizes are required, please call Comfortex and a customer service representative will assist you.
- * Most orders are processed, manufactured and shipped within 3-5 days.

Serving small facilities and large healthcare systems for more than 25 years.

 Comfortex

800-445-4007

*Being a Patient can be stressful enough.
Reduce Pressure ... Call Comfortex*

Comfortex Response Stretcher Cushion Order Form


Facility Name: _____
Contact Name: _____
City, State: _____
Department: _____
Phone: _____
E-mail: _____

Please fax completed form with purchase order to (507) 454-6581 or call Comfortex at (800) 445-4007 for assistance.

Note - Top Covers : Soft, shear reducing blue nylon is standard. Please indicate if LectoLite vinyl is preferred.


Comfortex will confirm your order and provide you with an estimated ship date via email....and thank you.

Standard Rounded Corners

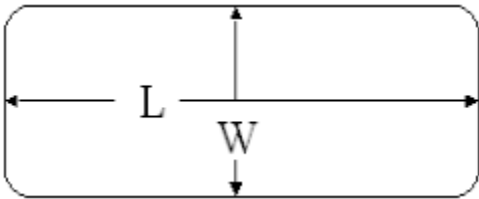


Length: _____
 Width: _____
 Quantity: _____
 Vinyl Preferred:

Standard Diagonal

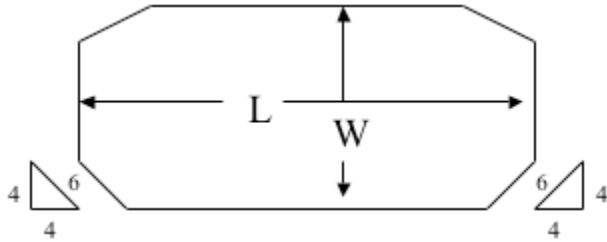


Length: _____
 Width: _____
 Quantity: _____
 Vinyl Preferred:

Head End Specific

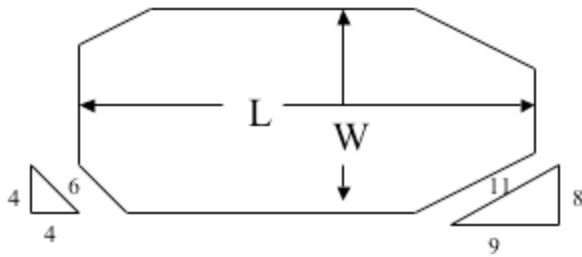

Length: _____
 Width: _____
 Quantity: _____
 Vinyl Preferred:

Custom (Add \$15.00 each)

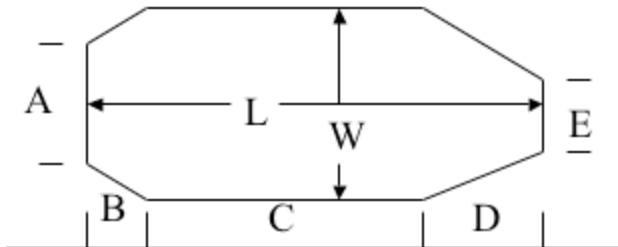

Length: _____
 Width: _____
 A: _____
 B: _____
 C: _____
 D: _____
 E: _____
 Quantity: _____
 Vinyl Preferred:

Response™

from  Comfortex

Stretcher Cushions

About Comfortex

We Work for You

"Comfortex recognizes that the success we experience will be in direct proportion to the benefits we provide our customers. Each design, each service and each action is guided by our mission to contribute to the quality of life for the patient, complement the efforts of the healthcare professional and enhance the financial strength of our Customers."

Andrea Terek, President

Craftsmanship

Materials like our Stratos covers and Engineered Polymer cores combine to make our surfaces comfortable years longer than any other. Technology blended with Craftsmanship ensures the surfaces you buy this year will perform like new long after. The patients you care for in the future deserve to be as comfortable as those cared for today.

Precision

From core fabrication to the stitches that secure the handles, Comfortex surfaces are crafted to last. The products that surround your patients are a reflection of Quality.

Comfortex Quality is a reflection of Excellence.

Customer Service

With more than 20 years in business and having placed hundreds of thousands of surfaces, Comfortex knows its best sales tool is a satisfied customer. Our representatives are committed to staying in touch to provide information, ensure benefits are experienced and that service and responsiveness exceed your expectations.

For quality & support...Comfortex reponds.

Hospitals Contact:

Comfortex®

Toll-Free:

(800) 445-4007

Comfortex,
an Avatar Enterprises,
Inc. company

1680 Wilkie Drive

Winona, MN 55987

Phone: (507) 454-6579

Fax: (507) 454-6581


E-Mail:

info@comfortexinc.com

Office Hours:

8 AM - 4:30 PM Mon-Thurs

8 AM - 2:00 PM Fri

 Comfortex

800-445-4007

*Being a Patient can be stressful enough.
Reduce Pressure ... Call Comfortex*