




A great idea...Reshaped

Design Features and Benefits of the Rest-Q VC

- 1.** Comfortex's Vapor Proof, Shear Reducing Top Cover is the Nation's most durable fabric in healthcare today. Other fabrics tend to fade, crack and can lead to core contamination.
- 2.** This layer of soft, high resiliency (HR) foam spans across the mattress to cushion boney prominences and ensure greater comfort, especially over the bed's exit points.
- 3.** Additional soft (HR) material further enhances comfort and pressure redistribution. This layer also incorporates extra supportive perimeter foam that overlays the bed's raised metal edging at the enter and exit points.
- 4.** Embedded within the core are supportive sections that underly the hips and shoulders. Between these support sections is a soft, pressure absorbing strip to help bridge the patient's pressure sensitive areas along the spinal column and sacrum. When side-lying, the trochanter rests bridged, over this center strip.
- 5.** A base layer of firm material helps ensure heavier patients will not bottom out. Both this firm base and the side-rail perimeters have relief cuts to allow the mattress to bend with the bed's position changes.
- 6.** Rather than slope the lower legs downward, which could cause swelling, the heels rest into the dual-density top foam and our ***descending diamond heel suspension system***. This heel suspension system absorbs the boney prominences of the heels while keeping the lower legs supported along the horizon of the mattress. The side-to-side chambers help accordion the foot end of the mattress if the bed frame is telescoped inward.

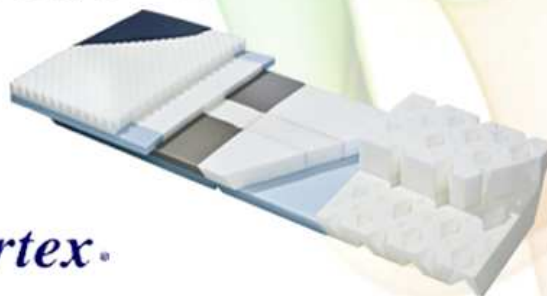

1.800.445.4007

www.comfortexinc.com

Wound Treatment...Pressure Sore Prevention, Support and Comfort


World class excellence supported by years of documented performance.

Pressure testing illustrates the Rest-Q's potential.


Multiple layers and zones of High Resiliency materials are crafted to fit all frames and clinical challenges.



Interface Pressure Imaging

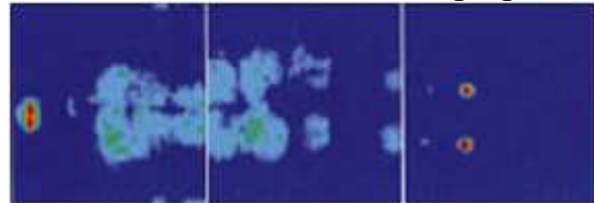

Low Airloss Mattress

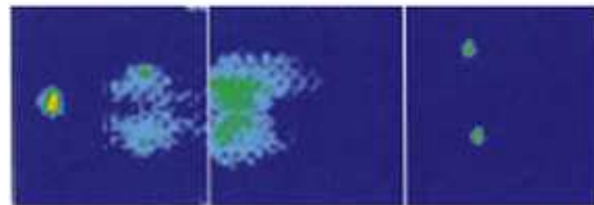

Rest-Q Mattress

Pain Reduction & Comfort


A German Institute of Health pain study compared the Comfortex Rest-Q® mattress and the Tempur® Memory Foam mattress to the standard mattress with 80 patients suffering chronic pain.

- Awoke as a result of pain: reduced by 52% on the Rest-Q, 21% on the Tempur mattress.
- Call lights reduced by 43% on the Rest-Q, 18 % on the Tempur mattress.
- Pain medications were 24% less on the Rest-Q and no reduction noted on the Tempur mattress.

Rated Most Comfortable

Standard			13	30	25	12
Tempur		2	18	7	8	6
Comfortex	5	13	9	6	5	1
	Very good	Good	Satisfactory	Adequate	Unsatisfactory	Poorly

Summary of Wound Change All Stages combined

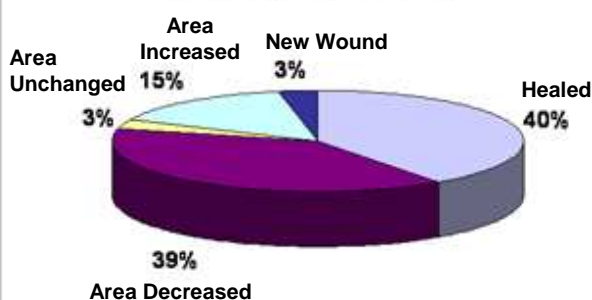

Treatment Capabilities (Independent Studies)

30 patients with pressure wounds were taken from more costly Group II surfaces and moved to the non-mechanical Rest-Q® mattress.

Nearly 80% of those wounds improved or healed.